Technical Data Sheet

InVivoPlus anti-mouse CD25 (IL-2Rα)



<u>Attention</u>: Use of this product constitutes an agreement to Bio X Cell's Terms and Conditions which are included with this product in print and can also be found at https://bioxcell.com/terms-and-conditions.

Lot Specific Information

Lot Number: Lot Specific* Volume: Lot Specific*

Concentration: Lot Specific* (generally 4 to 11 mg/ml) *

Total Protein: Lot Specific*

*This information will be noted on the certificate of analysis that ships with this product.

Product Information

 Catalog Number:
 BP0012

 Clone:
 PC-61.5.3

 Isotype:
 Rat IgG1, λ

Recommended Isotype Control(s): InVivoPlus rat IgG1 isotype control, anti-horseradish peroxidase

Recommended Dilution Buffer: InVivoPure pH 7.0 Dilution Buffer

Immunogen: IL-2-dependent cytolytic mouse T cell clone B6.1

Reported Applications: in vivo regulatory T cell depletion

Flow cytometry

Formulation: PBS, pH 7.0

Contains no stabilizers or preservatives

Endotoxin: <1EU/mg (<0.001EU/μg)

Determined by LAL gel clotting assay

Purity: >95%

Determined by SDS-PAGE

Sterility: 0.2 µm filtration

Production: Purified from cell culture supernatant in an animal-free facility

Purification: Protein A High Salt

Aggregation: <5%

Determined by SEC

RRID: AB_1107619
Molecular Weight: 150 kDa

Murine Pathogen Test Results

Mouse Norovirus: Negative, Mouse Parvovirus: Negative, Mouse Minute Virus: Negative, Mouse Hepatitis Virus: Negative, Reovirus Screen: Negative, Lymphocytic Choriomeningitis virus: Negative, Lactate Dehydrogenase-Elevating Virus: Negative, Mouse Rotavirus: Negative, Theiler's Murine Encephalomyelitis: Negative, Ectromelia/Mousepox Virus: Negative, Hantavirus: Negative, Polyoma Virus: Negative, Mouse Adenovirus: Negative, Sendai Virus: Negative, Mycoplasma Pulmonis: Negative, Pneumonia Virus of Mice: Negative, Mouse Cytomegalovirus: Negative, K Virus: Negative

Description

The PC- 6 1.5.3 monoclonal antibody reacts with mouse IL-2R α also known as CD25, Ly-43, p55, or Tac. IL-2R α is the 55 kDa ligand-binding subunit of the interleukin 2 receptor alpha chain. IL-2R α is expressed on activated mature T and B lymphocytes, thymocyte subsets, pre-B cells, and T regulatory cells. IL-2R α has been shown to play roles in lymphocyte differentiation, activation, and proliferation. Alone, the IL-2R α binds IL-2 with relatively low affinity however, when IL-2R α associates with IL-2R β (CD122) and the common gamma chain (CD132) the complex binds IL-2 with high affinity. The PC-61.5.3 antibody has been shown to inhibit the binding of IL-2 to both the low and high affinity IL-2 receptor forms. Additionally,

Bio X Cell, LLC Page 1 of 2

the PC-61.5.3 antibody is commonly used to deplete CD4+FoxP3+ T regulatory cells in vivo.

Storage

Store at the stock concentration at 4°C. **Do not freeze.**

It is not uncommon for a floccule or precipitate to appear during storage. The floccule is typically buffer salts precipitating out of solution or a small bit of protein aggregation. For information on how to remove floccules or precipitates see our FAQ's at https://bioxcell.com/fags.

Protocol Information

Since applications vary, each investigator should use the application references as a guide to help estimate the appropriate dose or concentration. The dose or concentration can be further optimized experimentally in a dose response or titration experiment.


Application References

For a complete list of references, visit https://bioxcell.com/bp0012?bxcs=9k1b3a#tab_references or scan the QR code below.



Binding Validation

Validation data shown below confirms that this clone binds to its target antigen. For lot specific binding validation data, e-mail technicalservice@bioxcell.com.

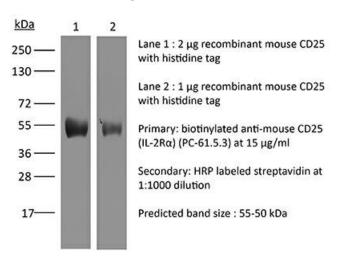

Bio X Cell, LLC

https://bioxcell.com +1-866-787-3444 customerservice@bioxcell.com Conditions: For research use only. Not for use in diagnostic or therapeutic procedures.

Not for resale.

Bio X Cell, Bio X Cell logo, and all other trademarks are the property of Bio X Cell, LLC © 2024 Bio X Cell, LLC

Bio X Cell, LLC Page 2 of 2